




1. Use the sizing gauge to select the correct mask size



2. Flip Head Strap forward and place mask on face. Pinch OmniClip™, slide up or down to position on patient's face

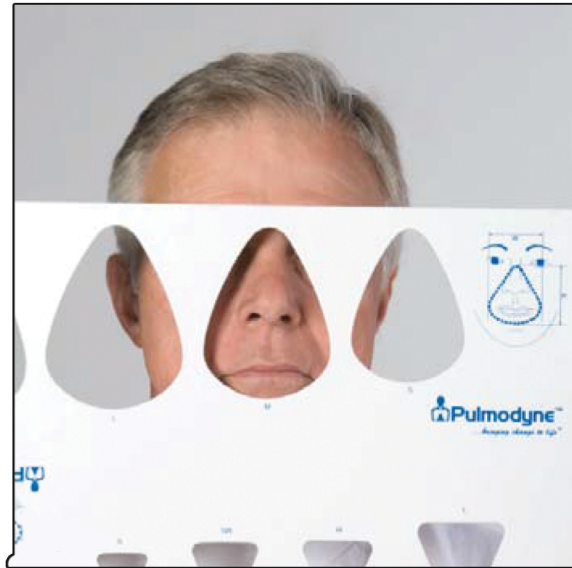

3. Flip Head Strap back over patient's head, bring tabs forward on "top" Head Strap and adjust equally to proper fit




4. Attach bottom two clips and repeat step #3



5. Connect circuit, turn on flow



6. Adjust OmniClip in and out, Head Strap and Mask for best fit. **Do Not Overtighten Head Strap**

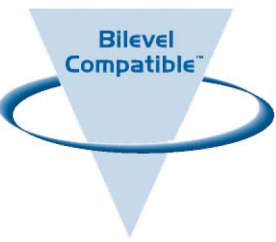

... bringing change to life™

Pulmonary and bringing change to life are registered trademarks of Pulmonary, Inc. BiTrac, OmniClip and Bilevel Compatible are trademarks of Pulmonary, Inc.

www.pulmonary.com